



Discovery DSC-MCAX (Mixed-Mode Cation Exchange) SPE Products


Base Silica: Irregular shaped, acid washed; 50µm particle size, 70Å pore diameter, 480m²/g specific surface area, 0.9cm³/g

Discovery DSC-MCAX SPE was developed, tested and quality controlled for pharmaceutical and clinical applications. The technology utilizes the dual mechanisms of both hydrophobic (reversed-phase) and electrostatic (ion-exchange) interaction to retain, basic, acidic, neutral and zwitterionic compounds from aqueous sample matrices. The packed bed contains both **octyl (C8)** and **benzene sulfonic acid (SCX)** bondings to offer one of the most selective SPE systems available.

Use Discovery DSC-MCAX SPE to:

- Achieve superior selectivity/sample clean-up when isolating basic/zwitterionic compounds.
- Experience greater and more reproducible recoveries of pharmaceutical compounds from biological samples
- Fractionate basic/zwitterionic compounds from acidic and neutral compounds.
- Adjust the protocol to permanently retain unwanted basic/ zwitterionic compounds while eluting acidic and neutral compounds of interest.

Recommended Generic Protocol:

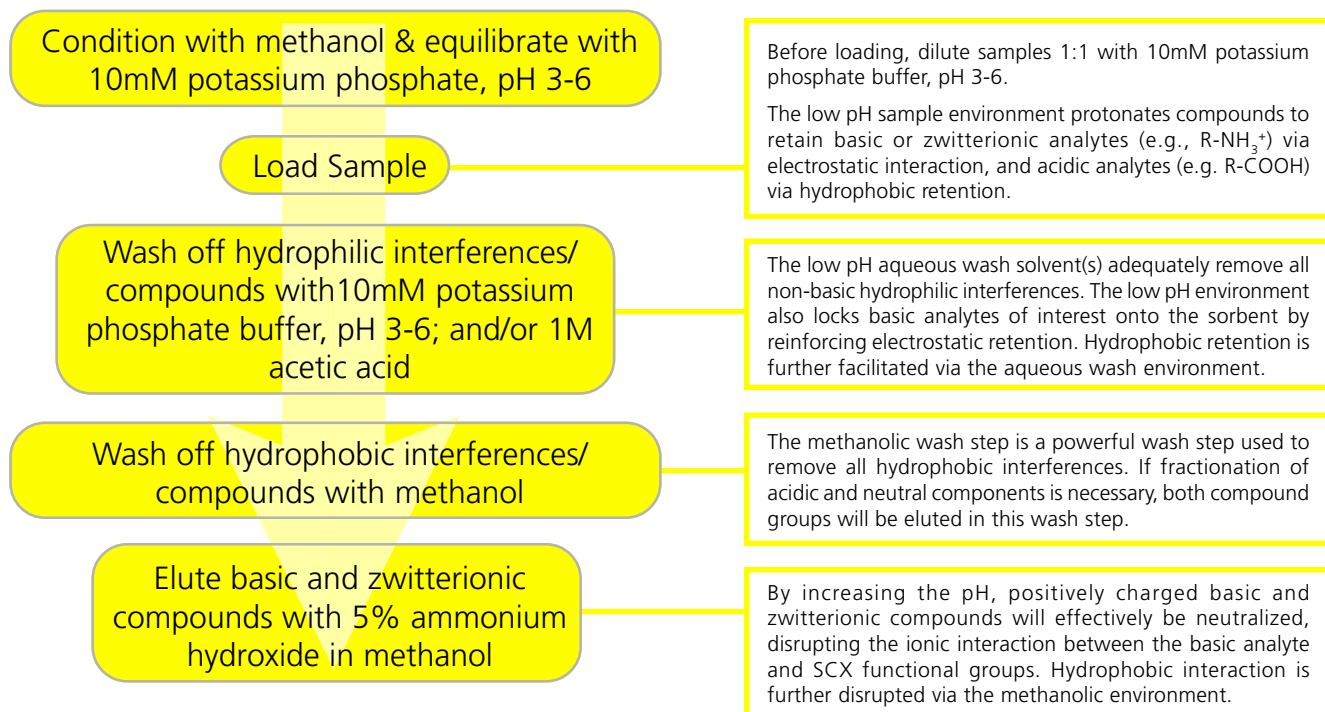


Here's the Proof:


Competitor Comparison of 2-Aminopyridine and Piroxicam from Human Urine

SPE Tube: Discovery DSC-MCAX, 100mg/3mL
 HPLC Column: Discovery HS F5, 15cm x 4.6mm, 5µm particle size
 Mobile Phase: 10mM potassium phosphate, pH 6:MeCN (85:15)
 Flow Rate: 2mL/min
 Temp.: 25°C
 Det.: UV, 220nm
 Inj. Vol.: 10µL


Spiked Urine Samples Before DSC-MCAX

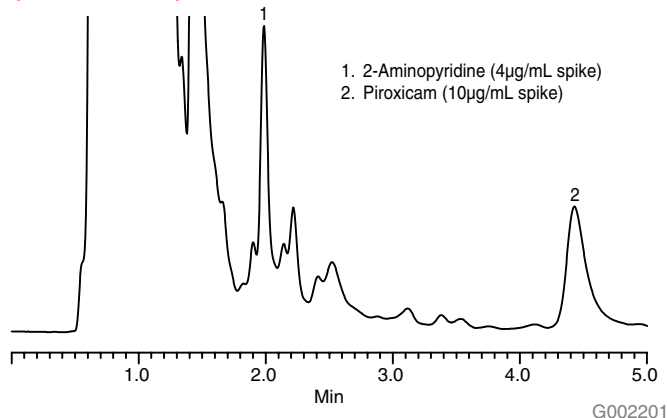


Spiked Urine Samples After DSC-MCAX

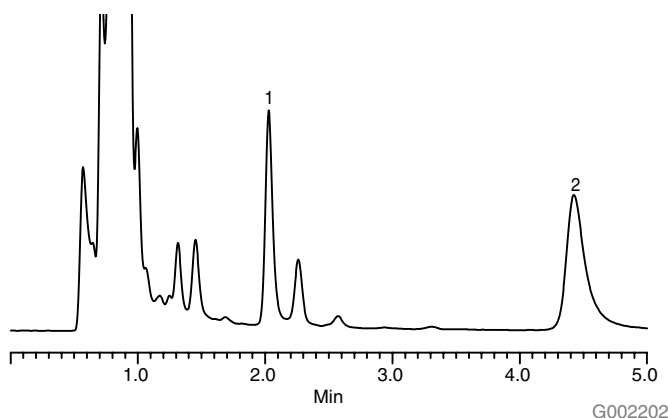

2-Aminopyridine

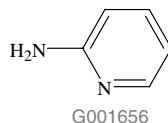

Piroxicam

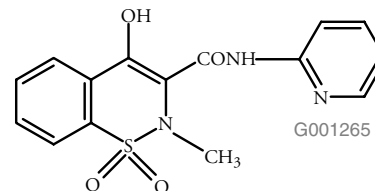

Efficiency of Absolute Recovery (n=4)

	% Recovery ± RSD	
	2-Aminopyridine	Piroxicam
Discovery DSC-MCAX	102 ± 3.5%	101 ± 1.2%
Leading Competitor A	30 ± 52.5%	98 ± 3.2%
Leading Competitor B	36 ± 24.2%	83 ± 4.3%

Ordering Information

Description	Qty./Pk	Cat. No.
Discovery DSC-MCAX		
50mg/1mL	108	52781-U
100mg/1mL	108	52782-U
100mg/3mL	54	52783-U
300mg/3mL	54	52784-U
300mg/6mL	30	52786-U
1g/6mL	30	52788-U
Discovery DSC-MCAX SPE 96-Well Plate		
25mg/well	1 ea.	575639-U
50mg/well	1 ea.	575640-U
100mg/well	1 ea.	575641-U

Trademark

Discovery is a trademark of Sigma-Aldrich Co.

sigma-aldrich.com/supelco

Order/Customer Service 800-247-6628, 800-325-3010 • Fax 800-325-5052 • E-mail supelco@sial.com
 Technical Service 800-359-3041, 814-359-3041 • Fax 800-359-3044, 814-359-5468 • E-mail techservice@sial.com
SUPELCO • 595 North Harrison Road, Bellefonte, PA 16823-0048 • 814-359-3441

ISO 9001
REGISTERED

We are committed to the success of our Customers, Employees and Shareholders through leadership in Life Science, High Technology and Service.

The SIGMA-ALDRICH Family SIGMA ALDRICH Fluka Riedel-deHaen SUPELCO

© 2003 Sigma-Aldrich Co. Printed in USA. Supelco brand products are sold through Sigma-Aldrich, Inc. Sigma-Aldrich, Inc. warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product(s) for their particular use. Additional terms and conditions may apply. Please see reverse side of the invoice or packing slip.



T403137
GEO